



**ARŞİMED MEDİKAL KOZM. SAN. TİC.
LTD. ŞTİ.**

Saray Köyü Organize Sanayi Bölgesi
Mevkii 4.Cadde No:5
YERKÖY/ YOZGAT - TURKIYE
Tel:+90 216 567 8 567 Fax:+90 216 567 8 367
info@arsimed.com.tr www.arsimed.com.tr

DC04-S PANAMERA DIALYSIS AND CHEMOTHERAPY CHAIR (5 Motors)



Arsimed Panamera Dialysis Chair is designed for dialysis, oncology, chemotherapy. The chair consists with 3 main parts. Back rest, seat part and leg rest. The chair makes back up&down, chair up&down, leg rest up&down, cramp support and trendelenburg&reverse trendelenburg movements by 5 DC motors with using hand control.

All adjustments done by hand control and nurse control unit. Cramp support is adjustable by motor for the cramps which happens during the dialysis.

Head part is adjustable for the patients' neck. Neck Cushion is given with the chair.

Chair is covered with cleanable, waterproof material and the chair is comfortable and orthopedic designed.

Cushion is suitable hardness and upholstery made by leatherette which is easy cleanable and not effected to skin. The chair is made from viscoelastic foam.

Viscoelastic foam, also known as memory foam, is a type of open cell, flexible polyurethane foam. It is body heat and weight sensitive material. It contours to the body shape, distributes the body pressure over a wide area and besides comfortable seating, allows for good blood circulation. The patient can also



ARŞİMED MEDİKAL KOZM. SAN. TİC. LTD. ŞTİ.

Saray Köyü Organize Sanayi Bölgesi
Mevkii 4.Cadde No:5
YERKÖY/ YOZGAT - TÜRKİYE
Tel:+90 216 567 8 567 Fax:+90 216 567 8 367
info@arsimed.com.tr www.arsimed.com.tr

benefit from the foam's capacity to redistribute weight and surface pressure, potentially reducing decubitus ulcers (bed sores) caused by aggravated pressure between the skin and bony areas of the body.

ABS covering designed to keep itself in the upholstery and it blocks the dust which goes into the upholstery. During the back rest movement, arm rests move with the back rest at the same time.

The arm rests are left-right rotatable and height adjustable. They are easily removable.

The chair can be produced with eating table mechanism.


1 pcs IV rod is standard.

Chair has a weighing system with 4 loadcell and 1 indicator. From indicator, patient weight can be seen. Weighing capacity is 300 kg.

NURSE CONTROL PANEL



HANDSET CONTROL PANEL

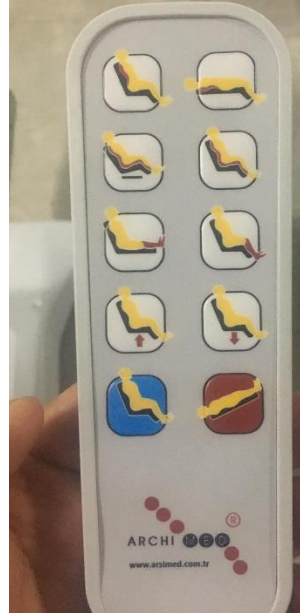



ARŞİMED MEDİKAL KOZM. SAN. TİC. LTD. ŞTİ.

Saray Köyü Organize Sanayi Bölgesi
Mevkii 4.Cadde No:5
YERKÖY/ YOZGAT - TÜRKİYE
Tel:+90 216 567 8 567 Fax:+90 216 567 8 367
info@arsimed.com.tr www.arsimed.com.tr

WEIGHING SCALE SPECS




Product Description

The T18 Weighing Indicator is designed to withstand heavy-duty working conditions thanks to its durable structure. It is resistant to shocks and vibrations, making it a strong choice for long-term use. Additionally, it is water and dust resistant, allowing it to be used outdoors or in challenging environments.

Technical Specifications

Product Code	T18
Accuracy Class	CLASS III
Display	1.6 LCD Display
Status Indicators	7 Status Indicator Signals
Keypad	Button Keypad
Resolution	0.6µV/e
Load Cell Input	-15mV ~ +15mV Input Range
Power Source	6V/4Ah Battery
Operating Temperature	-10°C ~ 40°C
ADC	A/D conversion with high stability and anti-interference ΔCZ
Supply Voltage	DC 5V
Load Cell Capacity	Load cell connection (4Wire) 1-4 units of 350 Ω or 8 units of 700 Ω
Calibration	With Reference Weight
Filtering	Digital Filtering
Weight Units	KG/LB
Serial Ports	Com1 – RS232
Body	Plastic Body
Weight	2.0 kg

Technical Drawing

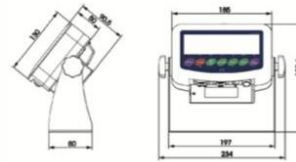

Application Areas

Protection Classification
Livestock Scales
Platform Scales
Truck Scales
Under
Over



**ARŞİMED MEDİKAL KOZM. SAN. TİC.
LTD. ŞTİ.**

Saray Köyü Organize Sanayi Bölgesi
Mevkii 4.Cadde No:5
YERKÖY/ YOZGAT - TÜRKİYE
Tel:+90 216 567 8 567 Fax:+90 216 567 8 367
info@arsimed.com.tr www.arsimed.com.tr

LEATHER COLORS





**ARŞİMED MEDİKAL KOZM. SAN. TİC.
LTD. ŞTİ.**

Saray Köyü Organize Sanayi Bölgesi
Mevkii 4.Cadde No:5
YERKÖY/ YOZGAT - TÜRKİYE
Tel:+90 216 567 8 567 Fax:+90 216 567 8 367
info@arsimed.com.tr www.arsimed.com.tr

PHOTOS FROM ONE OF OUR FINISHED PROJECTS

